


All-In-One Digital Microscope

The Digital Microscope has a built in 5" high resolution LED screen and offers a magnification range of 8x -25x. The included rechargeable battery allows you to operate the microscope without a power supply for hours thus

making it suitable for field projects. You can also connect this microscope to an external monitor via an HDMI output for viewing images on a larger screen. Images can be saved on the included Micro SD card. Can also be mounted onto a variety of stands to suit your application.

Features:

- Compact All-in-one system
- Save high resolution image/video files directly to microSD card
- Enhances productivity - eliminates eye strain and body fatigue
- 1080p Full HD Imaging at 60 FPS
- HDMI Output for additional monitor viewing
- software download included
- 1080p Full HD Output
- 5" built-in HD screen

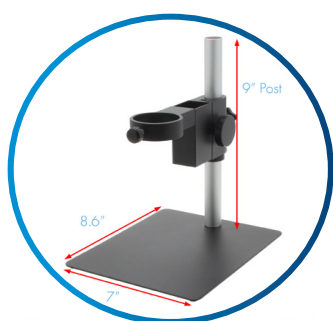

ezImage
Software
Included

Measure,
Share and
Annotate

8x-25x
Magnification

Adjustable
LEDs & Focus
Wheel


5" High
Resolution
LED Screen

1080P
Full HD

Micro SD
Card Slot

Process,
Save &
Share Images

Crystal Clear 5" Screen

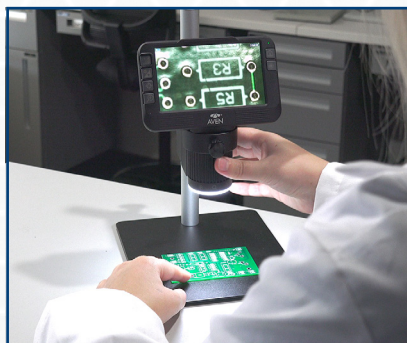

Software Download Included



Digital Microscope

Item #	#73-415
Magnification	8x-25x 40x-120x (on 24" Monitor)
Effective Pixels	1920(H) x 1080(V)
Image Sensor	2M CMOS Sensor
Working Dist.	45mm-95mm
Interface	Mini USB - Type A / HDMI
PC Software	EZImage X3
Illumination	Integrated Ring Light
Package Contents	Digital microscope, USB cable, Power supply, Snap shot remote control, 32gb MicroSD Card, 2600mAh 18650 L i-ion battery, User guide, Post Stand, Free software download